








Complete Blood Count (CBC) is a a routine blood test performed on pets to provide objective information about the general health status of an animal. The objective data obtained from a CBC can indicate the presence of many forms of diseases and can be helpful in monitoring ill animals undergoing therapy.

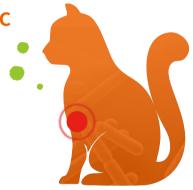

Complete Blood Count Applications



Providing evidence for diagnosis of inflammation or allergic diseases of animals, such as parasitosis, canine allergic bronchitis, cat asthma, allergic gastroenteritis, etc.



Aiding in diagnosing different kinds of anemia and polycythemia.



Aiding in diagnosing thrombocytosis. Primary thrombocytosis: essential thrombocythemia, myelodysplastic syndrome, etc.

Secondary thrombocytosis: high platelets caused by acute infection, hemorrhage, hemolysis, etc.

Species of Animal

Dog, cat, pig, cow, sheep, goat, horse, monkey, mouse, rabbit, pot-bellied pig, camel, buffalo, 2 self-programmed animals.





Fast Test

Utilizing double-channel test mode

No Residue

- No residue of hematology diluent or lysing reagent
- Cost effective

Humanization

- Sample volume of venous blood: 6 μL Less pain during blood collection
- Sample type: venous blood/capillary blood (pre-diluted mode), meeting different needs

Easy and Convenient Operation

- Large LCD touch screen
 High-definition color monitor, no need of external computer
- Built-in printer to fast print test results

Optimized Liquid System

Very low failure rate & excellent maintainability

Intelligent Floating Threshold

• Making the classification more accurate

BHA-3000 Vet

Automatic Hematology Analyzer

Technical Parameters

21 Reportable Parameters

WBC, NEU#, MXD#, LYM#, NEU%, MXD%, LYM%, RBC, HGB, MCV, RDW-CV, MCHC, MCH, RDW-SD, HCT, PLT, MPV, PDW, PCT, P-LCR, P-LCC

Sample Mode

Whole Blood (Venous/Capillary)
Pre-dilution

QC Methods

L-J QC, X-B QC

Screen Type

10.4 inch Resistive Screen

Weight

16 kg

Operating Conditions

18~30°C, ≤80% Relative Humidity, 75 kPa~106 kPa

Principles

Classic Electrical Impedance Method, Mature Photoelectric Colorimetry

Sample Volume

Venous blood: 6 μ L Capillary blood (Pre-dilution): 20 μ L

Data Storage

10000

Power Supply

100-240 V~ 50/60 Hz

Dimension

328 mm x 425 mm x 446 mm (L x W x H)



animal Add: No. 412 Wangxin Road, Luhe District, Nanjing, 211505, China Tel: +86-25-68568508/68568594

Tel: +86-25-68568508/685685 Fax: +86-25-57671919

E-mail: sales@getein.com.cn; overseas@getein.com.cn Web: www.getein.com









